bodynova

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Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Germany

PLEASE READ THE FOLLOWING SAFETY PRECAUTIONS CAREFULLY BEFORE USING.

1. IMPORTANT SAFETY INFORMATION

This symbol indicates a risk of electric shock inside the device that presents a risk of personal injury.

This symbol indicates important instructions for the user.

When using the Shiatsu Massage Cushion certain standard safety precautions, as with all other electrical devices, must be followed:

Use this device only inside, in a dry environment. Keep away from moisture and water. Do not use in the bath.

Keep the cushion and the cable away from heat sources.

The cable should not be bent, exposed or damaged in any way. Do not pierce the cable.

Disconnect the cable from the mains before cleaning, when not in use or when any problems arise.

Do not use the device for any other purpose than described in this manual.

Do not use while driving.

To avoid an electric shock keep the device, the cable and plug out of contact with water or any other liquids.

Do not use if the cushion, the cable or the plug are damaged in any way.

Do not use during thunderstorms.

If any problems arise while using the device, shut it off immediately and remove the plug from the mains.

Do not attempt to repair it yourself. Repairs to electrical devices should only be carried out by qualified personnel.

DISPOSAL:

Do not put the Shiatsu Massage Cushion 598BG in the regular waste disposal. It must be disposed of according to local regulations concerning electrical devices, see WEEE Guidelines (EU Guideline 2002/96/EG). For more information please contact your municipal waste removal company.

2. GENERAL USE


The handy Shiatsu Massage Cushion is a self-massager that provides an effective, deeply relaxing, kneading massage. The four rotating massage heads gently knead, loosen and soften the muscles along the spine. The optional warming function when switched on, delivers penetrating warmth deep into the tissue for a deeper relaxation.

The ergonomic shape of the cushion makes it perfect for massaging the lower back and along the spine. It is also excellent for using on the neck and shoulders, the legs and the feet. The small, compact form makes it an excellent travelling companion.

Simply rest the neck, feet or legs on the cushion for a relaxing massage or slide the cushion under the neck or back while sitting in the car, train or armchair.

The device is intended for personal use only and not for commercial use.

3. TECHNICAL DESCRIPTION

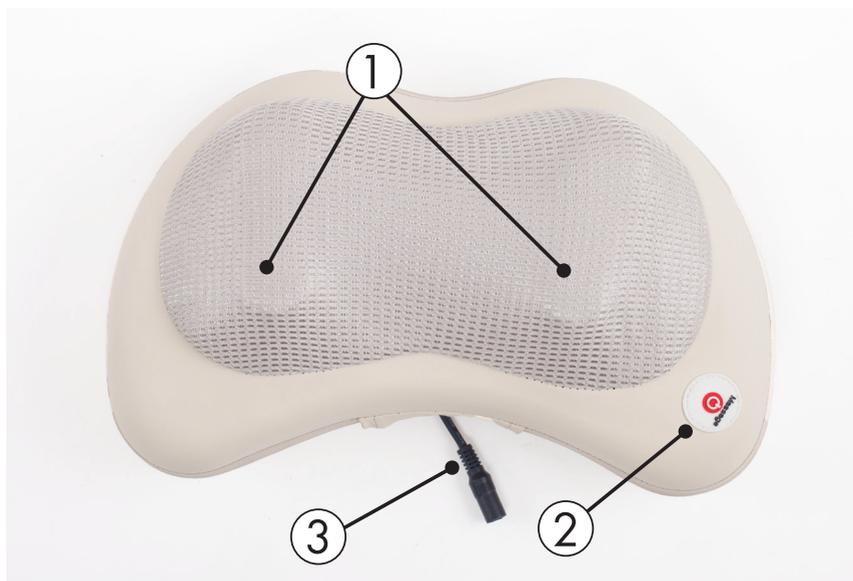1. Massage heads
2. on/off - control switch
3. Electrical connection

Included: 1 x 12V car adapter
 Power requirements:
 AC 230 V / DC 12V =2A

4. INSTRUCTIONS FOR USE

Before using the device, ensure that the cable is not damaged in any way. Connect the electrical cable to the connector on the device and connect the electrical plug to the mains.

While sitting comfortably on a chair, place the massage cushion between your lower back and the chair. Press the on/off switch once to turn the device on, it will begin massaging. Press the button again to activate the warming function. The optimal temperature of the warming function is reached after 6 minutes of use. Press and hold down the on/off switch to turn the device off. By increasing the pressure of the body on the massage heads you can vary the intensity of the massage effect.

The Shiatsu Massage Cushion 598BG is designed for use on the muscles along the spine. Simply place the cushion between the back of the chair and your spine. It can also be used on the muscles of the neck, the legs and the feet. Avoid using directly on bones or joints. Do not use longer than 15 minutes at a time. To avoid irritating the skin, be sure to use the device on top of your clothing and not directly on the skin. By discomfort or pain, stop using the massager immediately. If you have any health issues or are unsure, please consult your physician before use.

- The Shiatsu Massage Cushion 598BG has a built in automatic shut off and shuts down, after 15 min of use.
- After approximately 1 minute the rotating massage heads change direction.

5. CARE AND STORAGE

Before cleaning disconnect from the mains.

The cover of the massager is made of high-quality polyurethane and when needed, can be wiped clean with a damp cloth and a very small amount of a mild liquid detergent.

We also recommend cleaning the cover with DESTIX disinfection wipes, these are specially formulated for use on polyurethane.

Warning: aggressive detergents or disinfection fluids that contain a lot of alcohol can damage the material of the cover and should be avoided.

The massager should be stored in a dry room with temperatures between +10°C and +25°C.

6. CONTRAINDICATIONS AND PRECAUTIONS

The Shiatsu Massage Cushion 598BG is meant for use only as a self-massager and cannot replace proper medical treatment.

Use of the massager is contraindicated for all of the following:

- for all electrical implants (pacemakers, etc.) or metal implants
- for all types of injuries, open wounds and ruptured discs
- for pregnant women

Contact your physician if you have any questions or concerns with using the Shiatsu Cushion.

The device is only for consumer use and not meant for commercial or medical use.

Not to be used on pets or children and not for use by children.

The device is to be used only for the purposes described in this manual. The manufacturer accepts no liability for damages or injuries that occur as a result of improper use, or use contrary to the instructions contained in this manual.

7. WARRANTY

The manufacturer guarantees (for retail customers) this device to be free of defects in manufacturing for a period of 2 years from the date of purchase. Within this guarantee period, we will correct all defects due to material or manufacturing defects, either by means of repair or replacement, at our discretion.

Not included in this warranty are:

- damages caused by incorrect use
- normal wear and tear
- defects that do not impair the value or utility of the massager or do not impair these to a material degree

The warranty is void if repairs are attempted or performed by anyone not authorized by us. Any warranty claims within this 2-year period, whether it involves repairs or replacements, shall not extend the warranty beyond this period.

For warranty claims please send the device along with the purchase invoice or receipt to:
Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Germany

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Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Allemagne