bodynova

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Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Germany

PLEASE READ THE FOLLOWING SAFETY PRECAUTIONS CAREFULLY BEFORE USING.

1. IMPORTANT SAFETY INFORMATION

This symbol indicates a risk of electric shock inside the device that presents a risk of personal injury.

This symbol indicates important instructions for the user.

When using the Shiatsu Massagers certain standard safety precautions, as with all other electrical devices, must be followed:

Use this device only inside, in a dry environment. Keep away from moisture and water.
Do not use in the bath.

Keep the massager and the cable away from heat sources.

The cable should not be bent, exposed or damaged in any way.
Do not pierce the cable.

Disconnect the cable from the mains before cleaning, when not in use or when any problems arise.

Do not use the device for any other purpose than described in this manual.

Do not use while driving.

To avoid an electric shock keep the device, the cable and plug out of contact with water or any other liquids.

Do not use if the massager, the cable or the plug are damaged in any way.

Do not use during thunderstorms.

If any problems arise while using the device, shut it off immediately and remove the plug from the mains.

Do not attempt to repair yourself. Repairs to electrical devices should only be carried out by qualified personnel.

DISPOSAL:

Do not put the Shiatsu Massager 599BG in the regular waste disposal. It must be disposed of according to local regulations concerning electrical devices, see WEEE Guidelines (EU Guideline 2002/96/EG). For more information please contact your municipal waste removal company.



2. GENERAL USE

The Bodynova Shiatsu Massager is a self-massager that provides an effective, deeply relaxing, kneading massage. The massager is draped over the shoulders, like a shawl, allowing for a self-massage of the often tense, neck, shoulder and upper back area, which is otherwise very difficult to access yourself. It can also be used on the back, the legs and the feet. It can be used easily without any effort, simply resting the hands in the hand loops provided, allowing a deep and relaxing massage. In addition to the massaging action there is also a warming function that can be turned on and off, helping to ease and soften tension and stress.

With the hands or forearms in the hand loops provided, you can vary the depth and pressure of the massage by increasing or decreasing the amount of pressure applied through the hands.

The device is intended for personal use only and not for commercial use.

3. TECHNICAL DESCRIPTION

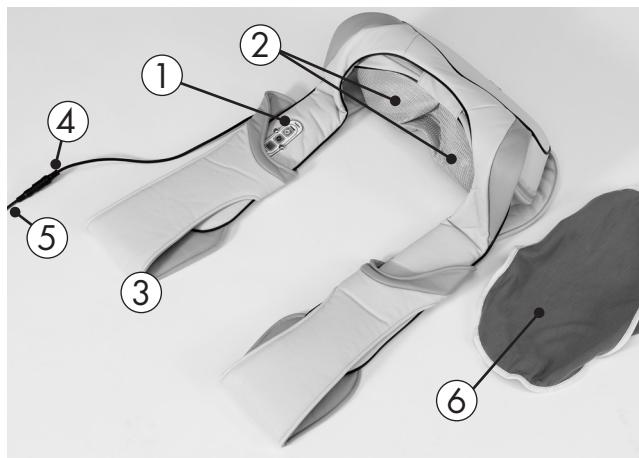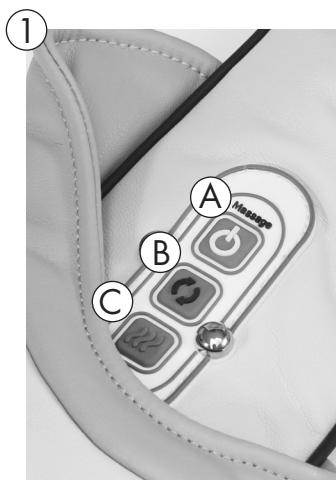1. Operating controls

- A: on/off switch
- B: switch to change the direction of rotation
- C: on/off switch for the warming function

2. Massage element

- 3. Hand loops
- 4. Connector for the AC adaptor
- 5. Electrical cable
- 6. Hygienic cover for massage element


Power requirements: AC 230 V / DC 12V =2A

4. INSTRUCTIONS FOR USE

Before using the device, ensure that the cable is not damaged in any way. Connect the electrical cable to the connector on the device and connect the electrical plug to the mains.

Drape the massager around your neck and over the shoulders, ensuring that it is centered. Turn the device on with the on/off switch. With your hands in the hand loops, you can adjust the device to find the optimal position and the pressure can be increased or decreased with the hands to achieve the ideal pressure for you. You can also turn the warming function on for an additional effect. By slightly adjusting the position of the massager you can reach different areas of muscular tension. The warming function can be switched on or off and the direction of the rotating heads can be changed. The optimal warming effect is reached after 6 minutes.

The Bodynova Shiatsu Massager 599BG is designed for use on the musculature of the neck, shoulders and upper back. The device can also be used for the muscles of the legs and feet, as well as the entire back. Do not use the device on bones or joints. Do not use longer than 15 minutes at a time. To avoid irritating the skin, be sure to use the device on top of your clothing and not directly on the skin. By discomfort or pain, stop using the massager immediately. If you have any health issues or are unsure, please consult your physician before use.

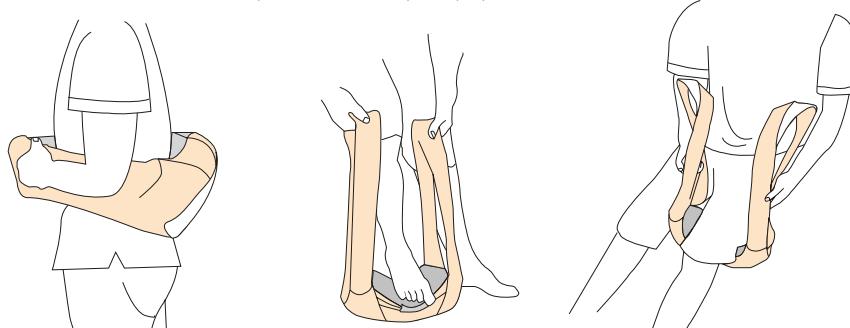Do not use longer than 15 minutes at a time. To avoid irritating the skin, be sure to use the device on top of your clothing and not directly on the skin. By discomfort or pain, stop using the massager immediately. If you have any health issues or are unsure, please consult your physician before use.

5. CARE AND STORAGE

Before cleaning disconnect from the mains. The cover of the massager is made of high-quality polyurethane and when needed, can be wiped clean with a damp cloth and a very small amount of a mild liquid detergent. We also recommend cleaning the cover with DESTIX disinfection wipes, these are specially formulated for use on polyurethane. Warning: aggressive detergents or disinfection fluids that contain a lot of alcohol can damage the material of the cover and should be avoided. The massager should be stored in a dry room with temperatures between +10°C and +25°C.

6. CONTRAINDICATIONS AND PRECAUTIONS

The Shiatsu Massager 599BG is meant for use only as a self-massager and cannot replace proper medical treatment.

Use of the massager is contraindicated for all of the following:

- for all electrical implants (pacemakers, etc.) or metal implants
- for all types of injuries, open wounds and ruptured discs
- for pregnant women

Contact your physician if you have any questions or concerns with using the Shiatsu Massager.

The device is only for consumer use and not meant for commercial or medical use.

Not to be used on pets or children and not for use by children.

The device is to be used only for the purposes described in this manual. The manufacturer accepts no liability for damages or injuries that occur as a result of improper use, or use contrary to the instructions contained in this manual.

7. WARRANTY

The manufacturer guarantees (for retail customers) this device to be free of defects in manufacturing for a period of 2 years from the date of purchase. Within this guarantee period, we will correct all defects due to material or manufacturing defects, either by means of repair or replacement, at our discretion.

Not included in this warranty are:

- damages caused by incorrect use
- normal wear and tear
- defects that do not impair the value or utility of the massager or do not impair these to a material degree

The warranty is void if repairs are attempted or performed by anyone not authorized by us. Any warranty claims within this 2-year period, whether it involves repairs or replacements, shall not extend the warranty beyond this period.

For warranty claims please send the device along with the purchase invoice or receipt to:
Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Germany

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Bodynova GmbH • Scheidtweiler Str. 19 • 50933 Cologne • Allemagne